


Heart Stimulating Electrodes



Mini ECG Electrode for Isolated Hearts

This special ECG electrode was designed to get the ECG (EG) from the surface of isolated hearts. The electrode can be used from the isolated mice heart up to the pig heart. In combination with the mini ball joint holders a precise fixation on the myocard is possible. Due to the flexible end these electrodes follow the beating heart. Connection: 4 mm banana connector, cable length 1 m.

Catalog No.	Model
CGS 8730.71	Mini ECG Electrode

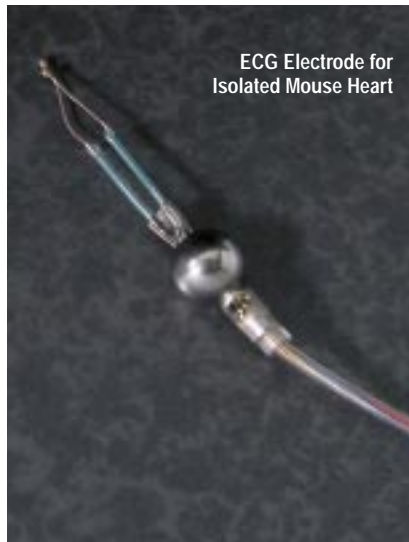

ECG Electrode for Isolated Heart

This miniature ECG electrode was designed to get the ECG from the surface of the isolated murine heart. It consists of two flexible silver wires with spheres fused to the end. In combination with the mini ball joint holders a precise fixation on the myocard is possible. Due to the flexible end these electrodes follow the beating heart. Connection: 4 mm banana connector, cable length 1 m.

Catalog No.	Model
CGS 8731.71 Descontinuado	ECG Electrode for Isolated Mouse Heart
CGS 8732.71	Connection Cable to Ph? ECGA Amp.



Coaxial Stimulation Electrode

This stimulation electrode is used for electrical stimulation of small heart or tissues. Through the coaxial construction the other electrode at zero potential acts as a screen so that stray fields are reduced and interference to recordings is limited (ECG, MAP). In combination with the mini ball joint holders a precise fixation on the myocard is possible.

Catalog No.	Model
CGS 8733.71	Coaxial Stimulation Electrode Mini
CGS 8734.71	Connection Cable to Stimulator

